





## All-day ECG Monitors

Hearts Deserve  
the Best Holters

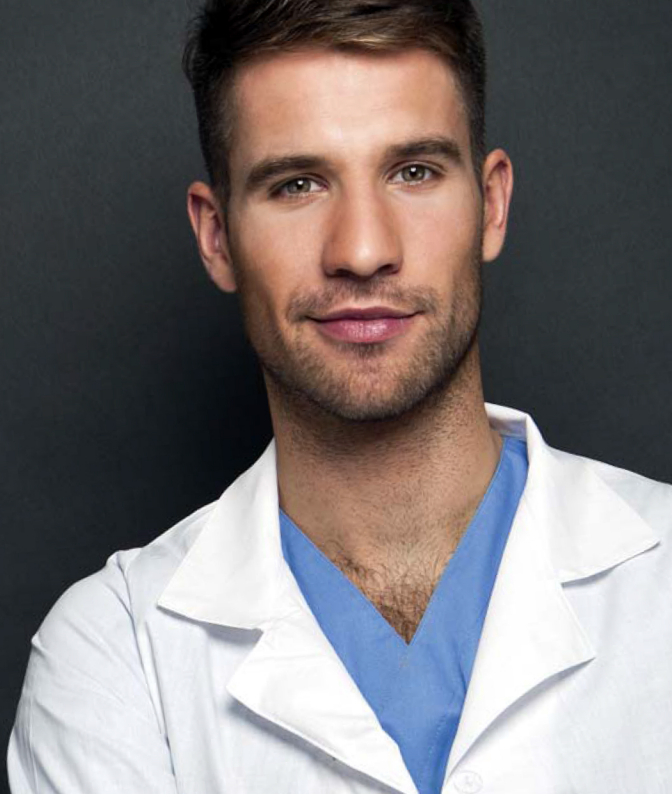

# CardiUP!12

3/12-channel Holter ECG

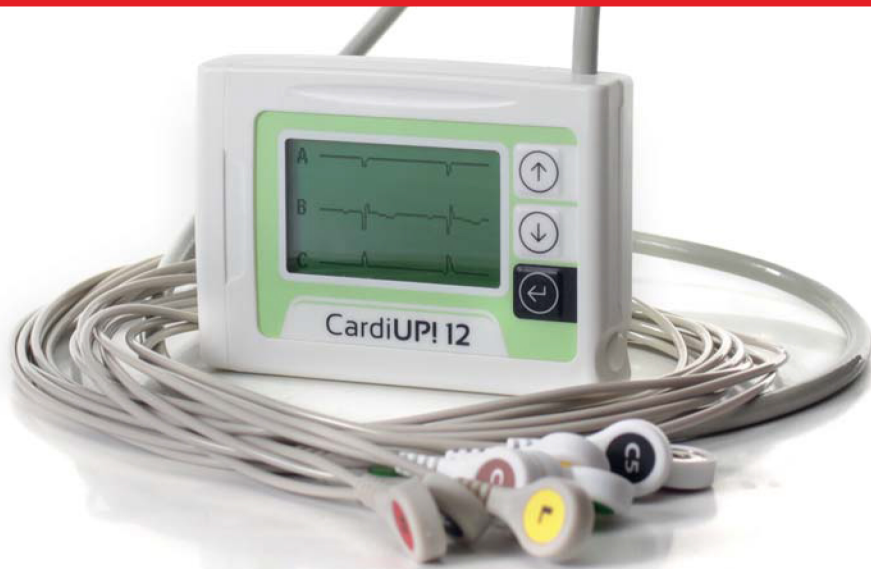

## CardiUP! benefits

- stay connected up to 7 days for better cardiovascular analysis
- 3/12 or 3 recording channels
- large LCD screen for real time ECG display
- event button
- both manual & software programming
- mini USB connection
- new, sensitive PM detection
- multifunctional carrying case
- free, all-inclusive software

## CardiUP!3

### 3-channel Holter ECG




# CardioMera

basic, multichannel Holter ECG

## CardioMera benefits

- 1-, 2-, 3- or 5-channel recording
- individual, easily replaceable patient cables
- small size & light weight
- event button
- free, all-inclusive software

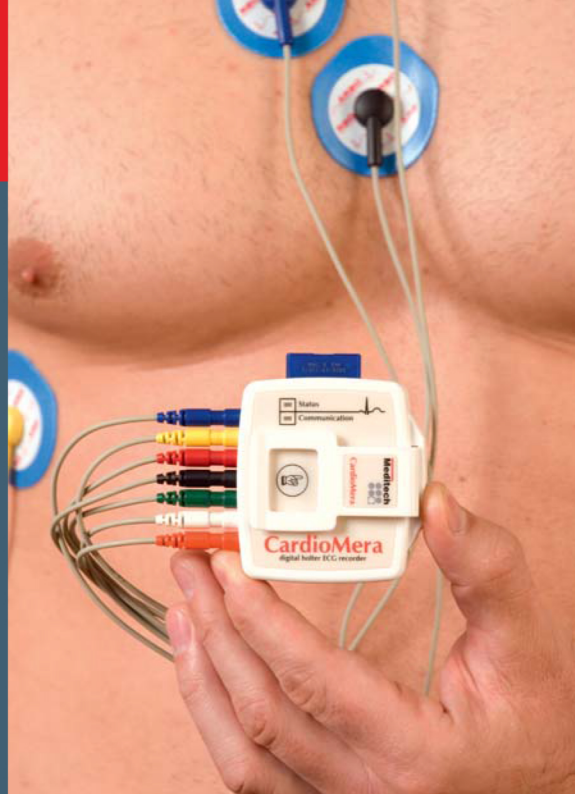

<b>Specifications</b>	<b>CardiUP!12</b>	<b>CardiUP!3</b>	<b>CardioMera</b>
dimensions	67*93*21 mm		66*59*17 mm
weight (excl. batteries)	92g		50g
memory	internal		external (a card reader unit is needed for retrieving data)
data transfer	USB 2.0 (mini USB)		USB optical cable
power supply	1 AA battery		1 AAA battery
LCD display	48*28 mm		-
number of keys	4 (event, up, down, select)		1 (event)
number of channels	3/12 independent channels (7 or 10 leads)	3 independent channels (7 leads)	1, 2, 3, 5 channels (7 leads)
operation	up to 168 hours of continuous ECG		up to 96 hours of continuous ECG
sampling, A/D resolution	300, 600 or 1200 Hz, 12 bit		600 or 1200 Hz, 12 bit
storage quality	300, 600 or 1200 Hz		150, 300 or 600 Hz
conformity certificates	CE		

## Holter ECG REPORT



### PATIENT DATA


name		Meditech	
ID	0001_0001	date of birth	2001.01.01.

### HOLTER SUMMARY


The measurement started	2018.03.09. 11:57	Total beats	96311	
Total measurement time	24 h 00 min	Normal	10553	(11,0%)
Valuable length:	24 h 00 min	VES	127	(0,1%)
		SVES	223	(0,2%)
Channels / sampling rate	12 / 300 Hz	PM	85408	(88,7%)
Total lead off / noise total	00 h 00 min	Other	0	(0,0%)
Pacemaker spike detection	enabled	Artifact	0	(0,0%)

16:46:48

A+  
HR=82

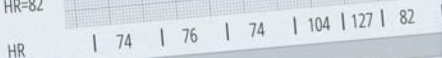

### HRV TACHOGRAM

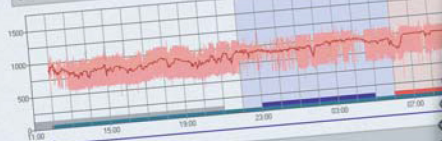

### HRV BAR CHART

SLEEP



## CardioVisions free software for ECG analysis

Specifications	CardiUP!12	CardiUP!3	CardioMera
accuracy	99,90%		
beat-by-beat ECG annotation	✓		
event distribution	✓		
ST analysis	✓		
QT analysis	✓		
PQ analysis	✓		
PM analysis	✓		
time domain HRV analysis	✓		
frequency domain HRV analysis	✓		
arrhythmia analysis	✓		
graphical and tabular display of data	✓		
customizable screens	✓		
customizable & printable report	✓		
personal & network use	✓		

## EXCLUSIVE DISTRIBUTOR IN AUSTRALIA

### Medusoft Pty. Ltd.

34 Devon Street, North Epping  
NSW 2121, Australia  
email: [info@medusoft.com.au](mailto:info@medusoft.com.au)  
[www.medusoft.com.au](http://www.medusoft.com.au)

## MANUFACTURER



Meditech Ltd • Mikszáth Kálmán u. 24. • 1184 Budapest, Hungary (EU)  
phone: +36 1 280 8232, +36 1 280 8233 • email: [meditech@meditech.eu](mailto:meditech@meditech.eu)  
[www.meditech.eu](http://www.meditech.eu)